

Regulatory information

CE Marking	Index of classification	LATEX
CE	I	

This stethoscope is made of plastic, for single-patient use only, in the prevention of nosocomial infections. Stethoscope with flat single head chestpiece made of ABS equipped with a non-chill ring and flexible ear-tips, insensitive to temperature changes. The binaural is made of polyamide.

1 – Device identification



2 – References


Reference	Designation
S20031	Stethoscope single head chestpiece Ø 34 mm Child
S20029	Stethoscope single head chestpiece Ø 45 mm Adult

3 – Use of the device - Operating instructions

- ⇒ Heart auscultation
- ⇒ Vascular auscultation
- ⇒ Pulmonary auscultation
- ⇒ Blood pressure monitoring

4 - Identification and characteristics of different subsets :

Identifier	Designation	Characteristic
A	Flat single head chestpiece	Aluminium
B	Screwed-type diaphragm	ABS
C	Ear tips	Flexible moulded PVC
D	Binaural with spring	Polyamide
E	Binaural	Polyamide
F	Tubing	PVC
G	Diaphragm	Terphane

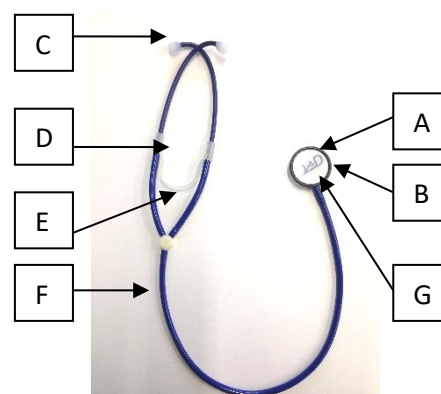

5 – Cleaning and disinfection

This stethoscope is not intended to be reused, it is not expected protocol for cleaning and disinfection.

6 - Packaging

Reference	Packaging unit
S20031	1
S20029	1
S20030	5

7 - Storage

Type of packaging	Storage area	Temperature	Humidity	Atmospheric pressure
Original packaging	Ventilated area	-10° to 40 ° c	30 to 40 %	500 to 1060 hpA



TECHNICAL CHARACTERISTICS

CT 022

Stethoscope STIMMER with single head chestpiece - Single-patient use

Index 6

8 - Warranty

This warranty provides assurance for the customer who purchases a D & D product that should the product fail to function to D & D published specifications during the term of this warranty, will either replace or repair. The product must be used in accordance with its labeling and may not be altered or subjected to misuse, abuse, accident or improper handling.

Dessillons & Dutrillaux

Zone Industrielle La Tuque
Tél : +33 5 53 48 30 66
technique@ddmedical.fr

47240 CASTELCULIER
Fax : +33 5 53 47 24 64
www.ddmedical.fr

Date of first EC marking certificate :
2004

Page 3 / 3