

The Best Choice of Accessories for Maximum Hygiene



Accessories Product range



The MEIKO solution for optimum allround-hygiene



Our motto „MEIKO – The clean solution“ encompasses more than just our portfolio of professional cleaning and disinfection appliances – it also includes the comprehensive range of MEIKO accessories. We believe in delivering optimum hygiene across the board, including surface and countertop materials in a selection of colours, machine accessories such as rack, holder and communication systems, plus a range of cupboards and shelf units.

Creating versatile and hygienic utility rooms



From single machines to fully-fitted utility rooms. Flexible, custom-made solutions including work surfaces, sinks, cupboards, shelf units and much more. MEIKO provides everything you need to make work easier and safer. High-quality stainless steel ensures optimum hygiene.

Easy to clean and available in a range of colours: High pressure laminate and CORIAN* panels



If you are looking for something special, we recommend high pressure laminates (HPLs). Manufactured from layers of paper impregnated with thermosetting resins, HPLs are both functional and attractive. Also available on request in other colours and designs at an additional cost.

A harmonious use of colour has been shown to have a positive effect on people's working environments. That's why we also offer our combined sluice units in a selection of colours on request.

Colour sample HPL

Melon yellow E2-31	Pastel green E23-21
Ivory E1-01	Sky blue E17-30
Sand E3-0	

Please contact us for full list of colours



CORIAN is a solid surface countertop material made from a blend of natural minerals and acrylic. Our combined sluice units can also be supplied in a colour of your choice. Wall and floor-mounted cupboard units can also be supplied in all RAL colours, spray painted in matte or glossy finish.

Colour sample CORIAN

Vanilla	Mojave
Raffia	Dusk

Please contact us for full list of colours

MEIKO machine accessories:
the right choice every time



The right choice every time: MEIKO machine accessories ensure you are prepared for every eventuality. From universal all-purpose racks to custom racks for special care utensils, we can help you find the right solution whatever your needs. Designed with users in mind, MEIKO accessories are practical and easy to use, providing reliable results and maximum hygiene.




Wash chamber door with lateral splash screens for use as slop sink
Cold water rinsing using the cleaning and disinfection machine's pump system.



SUH2 universal rack



Special holders for custom applications
Designed to pre-rinse out suction bottle and lid and dispose of the contents – all in a sealed environment.

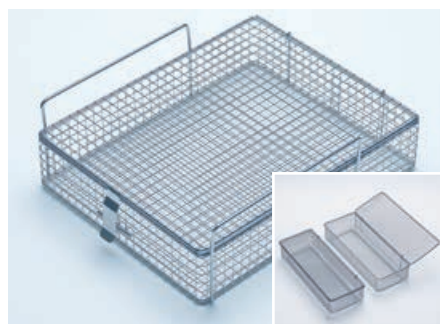


Basket for small items
With hinged lid. Designed to be inserted in the utensil holder (10 x 10 mm mesh size).

Basket insert
With hinged lid. Designed to be inserted in the basket for small items (4 x 4 mm mesh size).

MEIKO communication systems: Safeguarding hygiene standards



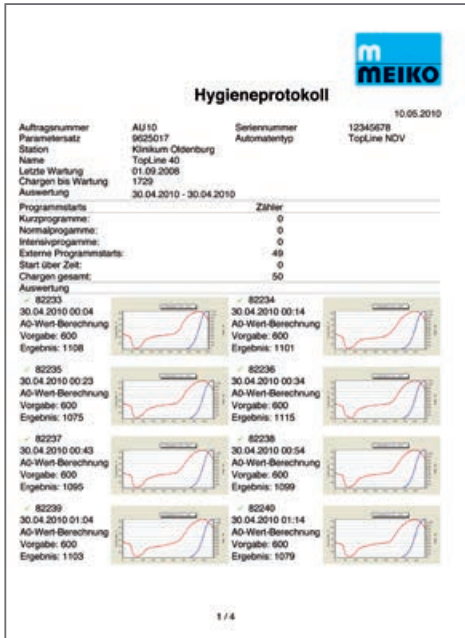

A seamless overview: *M-Commander InVision*

Stay up-to-date with comprehensive data on hygiene, logistical and technical aspects provided by MEIKO's *M-Commander InVision* analysis software. All the relevant information is transferred to your PC within seconds, so you can keep track of the machine's operating status clearly and easily. Compatible with MEIKO cleaning and disinfection appliances equipped with the *MIKE 2* electronic control system.

- A₀ value achieved
- A₀ value set
- date and time cycle started
- cycle duration
- machine identification no.
- cycle counter
- program no.
- temperature profile of disinfection phase
- online remote monitoring
- machine parameters
- and much more

Online disinfection monitoring with a wireless temperature data logger

Ask MEIKO to install a wireless temperature data logger in your cleaning and disinfection machine – and start monitoring your machine's disinfection performance online! Certified to DIN EN ISO 15883, this fast and reliable device is designed to ensure high standards of hygiene.

- Instantly browse process data on your laptop.
- The data is sent via LAN or radio transmission and checked by a MEIKO server within a matter of seconds.
- If everything is in order, a verification certificate is created automatically.
- This certificate is immediately sent to you for your files and can be used as evidence in the event of an inspection by the authorities.
- If one or more readings do not fall within the required parameters, the service engineer can immediately take the necessary steps on-site, thereby avoiding the need for another call-out.
- You retain permanent access to the appliance information and the measurement data.
- Obtain all the relevant data for your files within the scope of regular servicing

Clean solutions: From covers to feed pipes



Splash screen for slop sink

At rear and/or side. Optional rim flushing using cleaning and disinfection appliance or separate supply. Flushing rim on slop sink operated by touch sensor (dedicated supply or through appliance control system).



Cover with handle

Designed to prevent chemical odours. Available for soak sinks measuring 500 mm x 400 mm x 300 mm.



Single-lever mixer tap

Chrome plated, with swivel spout, easy to operate.



Lever waste valve

For sinks from size 500 mm x 400 mm x 300 mm, with optional fill level mark.



Cleaning agents and detergents: DOYEN SK 22 E special water softener and cleaning agent (also available with apple scent)

Prevents lime-scale deposits and rinses washware without leaving drip marks. Guarantees long-lasting cleaning and disinfection results. (5 kg container)

DOYEN R 100 bedpan cleaning agent

Mildly alkaline, powerful yet gentle cleaning action to combat heavy soiling. (5 kg container)

DOYEN RS 120 bedpan cleaning agent

Acidic, reliably breaks up fat and ointment residues, can be used at all levels of water hardness. (5 kg container)



Central drain piping

Central connection point under slop sink.



Feed pipes for mixer taps

Central connection for cold water and hot water (includes plumbing wall).



Stainless steel countertop with various different designs of wall edging

e.g. with straight upstand along wall or with hygiene radius (R15) and sloping upstand along wall.

Cupboards and shelves: Keeping utility rooms neat and tidy



Cupboards

e.g. cupboard unit with hand wash basin, 2 doors and 4 drawers.
Available in a range of sizes.



Regal TUS

For disinfected urine bottles, bedpans and receptacles, with storage rack or drip tray. Available in the following widths: 500 mm and 1,400 mm.

Wall shelving

Rack with height-adjustable shelves available in all shelf lengths. Shelf depth: 350 mm.



Cupboards without doors


Cupboard units with one intermediate shelf for quick and easy storage. Available in the following widths: 500 mm, 1,000 mm and 1,200 mm.

Cupboards with doors

Cupboard units with hinged doors and one intermediate shelf. Available in the following widths: 500 mm, 600 mm, 800 mm, 1,000 mm, 1,200 mm and 1,500 mm. (Optional: light screen with lamp)

Rack

For vertical storage of bedpans in wall-mounted cupboard units.

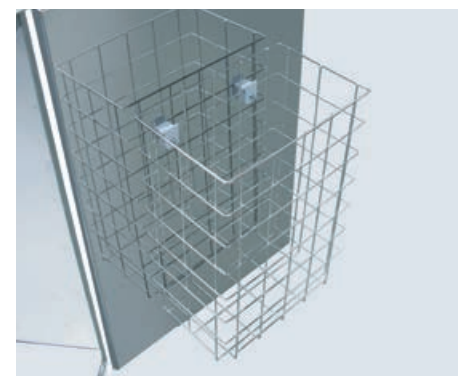

Lockable cupboard door

Prevents unauthorised access.

Rubbish bin

Stainless steel, with lid. Incorporated in lower cupboard unit for convenient disposal of small items of waste.


Holder for rubbish bag (26 liters)

Incorporated in lower cupboard unit for convenient disposal of small items of waste.



The clean solution

Manufacturer:
MEIKO Maschinenbau GmbH & Co. KG · Englerstr. 3 · 77652 Offenburg, Germany
Phone +49 781 203-0 · Fax +49 781 203-1179
www.meiko-global.com · info@meiko-global.com

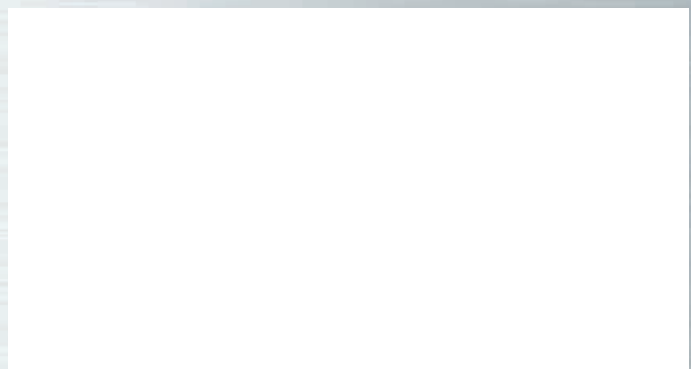

740.091.0X.XX.XX/EN/

We reserve the right to amend specifications as part of our product improvement process.